


Journal Pre-proof

Lack of sensitivity of an IVD/CE-labeled kit targeting the S gene for detection of SARS-CoV-2

Lucie-Marie Matzkies, Eva Leitner, Evelyn Stelzl, Karoline Assig, Michael Bozic, David Siebenhofer, Maria E. Mustafa, Ivo Steinmetz, Harald H. Kessler

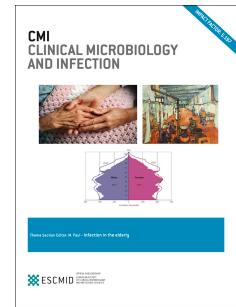

PII: S1198-743X(20)30382-7

DOI: <https://doi.org/10.1016/j.cmi.2020.06.036>

Reference: CMI 2125

To appear in: *Clinical Microbiology and Infection*

Received Date: 15 May 2020

Revised Date: 29 June 2020

Accepted Date: 30 June 2020

Please cite this article as: Matzkies L-M, Leitner E, Stelzl E, Assig K, Bozic M, Siebenhofer D, Mustafa ME, Steinmetz I, Kessler HH, Lack of sensitivity of an IVD/CE-labeled kit targeting the S gene for detection of SARS-CoV-2, *Clinical Microbiology and Infection*, <https://doi.org/10.1016/j.cmi.2020.06.036>.

This is a PDF file of an article that has undergone enhancements after acceptance, such as the addition of a cover page and metadata, and formatting for readability, but it is not yet the definitive version of record. This version will undergo additional copyediting, typesetting and review before it is published in its final form, but we are providing this version to give early visibility of the article. Please note that, during the production process, errors may be discovered which could affect the content, and all legal disclaimers that apply to the journal pertain.

© 2020 Published by Elsevier Ltd on behalf of European Society of Clinical Microbiology and Infectious Diseases.

Research note

Lack of sensitivity of an IVD/CE-labeled kit targeting the S gene for detection of SARS-CoV-2

Lucie-Marie Matzkies¹, Eva Leitner¹, Evelyn Stelzl¹, Karoline Assig¹, Michael Bozic¹, David Siebenhofer¹, Maria E. Mustafa², Ivo Steinmetz¹, Harald H. Kessler¹

¹Diagnostic & Research Institute of Hygiene, Microbiology and Environmental Medicine, Medical University of Graz, Graz, Austria

²Medical - chemical laboratory Dr. Mustafa, Dr. Richter OG, Salzburg, Austria

Correspondence: Ivo Steinmetz (ivo.steinmetz@medunigraz.at)

Neue Stiftingtalstraße 6, 8010 Graz, Austria

Telephone: 0043316-385-73700, Fax: 0043316-385-79647

Keywords: COVID-19, VIASURE, cobas, real time RT-PCR, diagnostics

Abstract

Objectives: New molecular tests for SARS-CoV-2 are rapidly launched in response to the COVID-19 pandemic. The aim of this study was to evaluate the analytical and the clinical performance of the VIASURE SARS-CoV-2 S gene RT-PCR Kit on the BD Max™ system and to compare results with those obtained with the cobas® SARS-CoV-2 test on the cobas® 6800 system.

Methods: For testing the analytical performance, reference material was used. Clinical samples (n=101) obtained from patients with symptoms compatible to COVID-19 were studied. Oro- and nasopharyngeal swabs were collected by using either ESwab™ or UTM™ collection systems.

Results: When the analytical performance was evaluated, the sample containing the lowest SARS-CoV-2 concentration tested negative with the VIASURE test while results obtained with the cobas® test were found to be concordant with the results expected. Six out of the 101 clinical samples (5.9%) showed an inhibition with the VIASURE test. When analyzing the remaining 95 clinical samples, 27 were found to be negative with both assays. Of 68 samples positive with the cobas® test, the VIASURE test missed 21 (30.9 %) samples. All of those 21 samples had shown Ct values ≥ 31 with the cobas® 6800 system. None of the samples tested positive with the VIASURE test and negative with the cobas® test.

Conclusions: The VIASURE test was impaired by a lack of sensitivity and a relatively high number of invalid results. When using the VIASURE test for routine testing, a significant number of COVID-19 positive samples would have been missed.

Introduction

The COVID-19 pandemic causes an extremely high demand for molecular SARS-CoV-2 testing of respiratory samples. This situation puts high pressure on manufacturers of molecular diagnostics to place new assays rapidly on the market to meet this expectation. Both laboratories in Graz and Salzburg involved in this study, are major COVID-19 diagnostic institutions in the Federal State of Styria and Salzburg in Austria, respectively. We have been testing for SARS-CoV-2 since the end of January 2020 and initially used the real time RT-PCR protocol established by Corman et al. [1]. Due to the significant increase of sample numbers, we introduced commercially available assays including the fully automated IVD/CE-labeled (and FDA-approved) cobas® SARS-CoV-2 test on the cobas® 6800 system (Roche Molecular Diagnostics, Branchburg, NJ). The analytical performance of this assay in clinical specimens was recently evaluated[2]. The assay utilizes a two-target RT-qPCR system targeting the ORF1 (target 1), a non-structural region that is unique to SARS-CoV-2, and a conserved region in the structural protein envelope E gene for pan-164 Sarbecovirus (target 2) detection [2]. In March 2020, the IVD/CE-marked VIASURE SARS-CoV-2 S gene RT-PCR Kit from CerTest Biotec (Zaragoza, Spain) targeting the S gene was launched [3]. This test has been designed for the BD MAX™ diagnostic platform (Becton Dickinson, Franklin Lakes, NJ) also integrating nucleic acid extraction, purification, amplification and detection. The latter platform has been widely used for the detection of bacterial and viral pathogens from clinical samples [4, 5] , including our diagnostic laboratories. As part of ongoing efforts to optimize the workflow of our COVID-19 diagnostics with different assays, we aimed to introduce the VIASURE SARS-CoV-2 S gene RT-PCR to be run on the BD MAX™ platform. Here we report on the results of independent

validations of this newly CE-marked automated SARS-CoV-2 S gene RT-PCR test from CerTest Biotec in our two laboratories.

Methods

In order to evaluate the clinical performance of the VIASURE test at the laboratories in Graz and Salzburg, cohorts of cobas® 6800 positive clinical samples with Ct values ranging from 20.82 to 37.85 and negative samples were compiled and retested by the VIASURE test according to the manufacturers' instructions. Oro- and nasopharyngeal swabs were either received from hospitalized patients, health care workers, or symptomatic individuals in the community presenting symptoms compatible with COVID-19. All samples were collected by using either ESwab™ or UTM™ collection systems (COPAN spa, Brescia, Italy). In total, 101 cobas® positive and negative samples were additionally analyzed with the VIASURE test. For this, 200 µl of sample transport medium were retested within 48h (stored at 2-8°C) after sampling by using VIASURE test lot number NCO124-004 and NCO124-010 in Graz and Salzburg, respectively.

Samples of the Quality Control for Molecular Diagnostics (QCMD) 2020 Coronavirus Outbreak Preparedness EQA Pilot Scheme (<https://www.qcmd.org/>) were tested with both the VIASURE SARS-CoV-2 S gene RT-PCR test (lot no. NCO124-004) on the BD MAX™ platform and the cobas® 6800 SARS-CoV-2 RT-PCR assay (lot no.G06568) on the cobas® 6800 platform according to the manufacturers' instructions at the Graz laboratory. The panel included samples containing different concentrations of SARS-CoV-2, coronavirus NL63, and coronavirus OC43 in transport medium as well as blank transport medium.

The study was carried out in accordance with The Code of Ethics of the World Medical Association (Declaration of Helsinki). Only pseudonymized left over specimens from routine

diagnostics were used to validate the performance of the VIASURE SARS-CoV-2 *S gene* Real Time PCR Detection KIT with no informed consent sought.

Results

Performance of the VIASURE SARS-CoV-2 *S gene* RT PCR Kit with samples of an external quality assurance (EQA) pilot scheme

Results obtained with the cobas® 6800 SARS-CoV-2 test were consistent with the results expected. With the VIASURE SARS-CoV-2 *S gene* RT PCR Kit, all samples containing SARS-CoV-2 concentrations between 2,000 and 200,000 copies/ml tested positive; however, the sample containing the lowest concentration (200 copies/ml) could not be detected.

Performance of the VIASURE SARS-CoV-2 *S gene* RT PCR Kit with clinical samples

Six out of the 101 samples (5,9%) showed an inhibition (unresolved result) with the VIASURE test (5 cobas®-positive and one cobas®-negative sample) and were excluded from further analysis. When analyzing the remaining 95 clinical samples, 27 were found to be negative with both assays. From 68 samples positive with the cobas® SARS-CoV-2 test, the VIASURE test missed 21 (30.9 %) samples (Figure 1). Samples found to be negative with the VIASURE kit but positive with the cobas® kit showed significantly higher cobas® 6800 cycle threshold (Ct) values than those testing positive with both of the kits ($p < 0,05$, t-test). The VIASURE test missed 52,5% (21/40) of samples positive with the cobas® test showing a Ct value of ≥ 31 . All samples that tested negative with the cobas® test were also found to be negative with the VIASURE SARS-CoV-2 test.

The false negative results with VIASURE test from both laboratories were communicated independently to BD. As proposed by the manufacturer, the Salzburg laboratory investigated whether sensitivity may be improved with a higher sample input volume and at the same time an adjusted higher volume for extraction according to manufacturer's instructions. Increasing the input volume with 21 samples from 200 µl to either 500 µl or 750 µl resulted in 4 (2 negative and 2 inhibited) and 2 (1 negative and 1 inhibited) discrepant test results, respectively, compared to the cobas® test.

Pharyngeal swabs with cobas® 6800 SARS-CoV-2 Ct values ≥ 31 are common in COVID-19 patients

We then aimed to evaluate the potential impact of false-negative results obtained in our compiled COVID-19 samples with the VIASURE assay (Figure 1) for routine COVID-19 testing. All Ct values of throat swabs that had tested positive for SARS-CoV-2 RNA with the cobas® test at the Graz lab between March 27th and April 2nd 2020 were analyzed. We found that 50.9 % (298/586) of positive samples showed Ct values ≥ 31 (Figure 2). Those samples included 224 newly diagnosed patients. Among those were 162 newly diagnosed patients in the community and 62 newly diagnosed patients from hospitals. 74 samples were derived from previously diagnosed patients. When taking into account, that 52.5% (21/40) of cobas® test positives with Ct values ≥ 31 were not detected by the VIASURE assay (Figure 1), one could estimate that approximately 26.7% of our COVID-19 samples would have been missed by this assay.

Discussion

In contrast to SARS-CoV-1, the causative agent of COVID-19, SARS-CoV-2, replicates efficiently in the upper respiratory tract leading to high viral loads in throat swabs during the first week of symptoms [6]. However, in most patients with mild disease viral load rapidly declines during the second week of illness [6]. Therefore, the high proportion of throat swabs with a relatively low viral load (Ct values ≥ 31) derived from hospitalized patients and patients in the community in our study is not surprising. COVID-19 respiratory tract samples with Ct values ≥ 31 were also observed by others in different SARS-CoV-2 RT PCR assays [7, 8]. Apart from the dynamics of the viral load in a short period of time the sampling procedure of oropharyngeal or nasopharyngeal swabs is likely to add to highly variable viral loads detected in those specimens. Therefore, highly sensitive real time RT-PCRs are most important for any case detection and preventive measures. For COVID-19 RT-PCRs intended for use by health professionals, the Directive 98/79/EC devices (IVD) is currently applied [9]. According to this Directive, a CE-mark can be affixed to the test, when the manufacturer declares that the requirements of the Directive are fulfilled (declaration of conformity) without involving a notified body. The CE-marked VIASURE SARS-CoV-2 S gene Real Time PCR Kit for use with the BD MAXTM is an automated system in which cartridges are used for the full processing of single samples [3]. According to the Instructions for Use, the clinical sensitivity and specificity of the VIASURE SARS-CoV-2 S gene RT-PCR Kit was tested using 4 nucleic acids isolated from SARS-CoV-2 positive respiratory samples and 15 respiratory samples from patients with clinical suspicion of COVID-19 disease or other respiratory diseases. From this panel 4 SARS-CoV-2 positive samples were detected and results were confirmed by two manual assays [3]. Direct comparison of all test results obtained with this limited number of specimens to a reference method was not reported [3].

Our study revealed a high percentage of false negative Covid-19 samples with the VIASURE SARS-CoV-2 S gene Real Time PCR Kit, although CE marked [9] when validated with an appropriate target population and reference method. The small sample size is a limitation of our study, but is also due to an extreme shortage of SARS-CoV-2 assay reagents during the study period. Although the VIASURE Real Time PCR Cartridges are provided with an internal PCR inhibition control, positive controls with the SARS-CoV-2 target gene are not provided. A possible reason for missing detection may be sensitivity issues related to the S gene target which might be less suitable for SARS-CoV-2 diagnostics or nucleic acids extraction and removal of inhibitors. The latter two being essential for reliability of results in routine molecular diagnostics [10]. Sensitivity issues are also of major importance when it comes to pooling of samples for SARS-CoV-2 testing. As reported recently, pooling affects the sensitivity leading to an increase of the Ct value and the false-negative rate [11].

Given the high numbers of newly CE-marked COVID-19 RT-PCR assays, the European Commission very recently communicated guidelines on COVID-19 in vitro diagnostic tests and their performance [12]. Although the Directive 98/79/EC defines a number of test evaluation data which have to be provided by the manufacturer to document test performance, the European Commission recommends to carry out additional clinical validation of COVID-19 diagnostic tests by comparison with a reference method including scientific peer-review [12]. The results of our study strongly support this approach before introducing a COVID-19 test into clinical routine.

Authors' contribution

LM: conceptualized the laboratory work, conducted laboratory tests, analyzed data, prepared figures and wrote the manuscript; EL: conceptualized the laboratory work, analyzed data; ES: analyzed data; KA: analyzed data and prepared figures; DS: conducted laboratory tests; MB: conducted laboratory tests; MM: conceptualized the laboratory work, supervised laboratory tests and analyzed data; IS: conceptualized the overall study, analyzed data and wrote the first manuscript draft; HK: analyzed data and edited the manuscript

1. Corman VM, Landt O, Kaiser M, Molenkamp R, Meijer A, Chu DK, et al. Detection of 2019 novel coronavirus (2019-nCoV) by real-time RT-PCR. *Euro Surveill* 2020;25.
2. Poljak M, Korva M, Knap Gašper N, Fujs Komloš K, Sagadin M, Uršič T, et al. Clinical Evaluation of the cobas SARS-CoV-2 Test and a Diagnostic Platform Switch during 48 Hours in the Midst of the COVID-19 Pandemic. *J Clin Microbiol* 2020;58.
3. Instructions. 444212 BR. VIASURE SARS-CoV-2 S-gene Real Time PCR Detection Kit. (March 2020).
4. Knabl L, Grutsch I, Orth-Holler D. Comparison of the BD MAX(R) Enteric Bacterial Panel assay with conventional diagnostic procedures in diarrheal stool samples. *Eur J Clin Microbiol Infect Dis* 2016;35:131-6.
5. Grimmond T, Neelakanta A, Miller B, Saiyed A, Gill P, Cadnum J, et al. A microbiological study to investigate the carriage and transmission-potential of *Clostridium difficile* spores on single-use and reusable sharps containers. *Am J Infect Control* 2018;46:1154-9.
6. Wölfel R, Corman VM, Guggemos W, Seilmaier M, Zange S, Müller MA, et al. Virological assessment of hospitalized patients with COVID-2019. *Nature* 2020;581:465-9.

7. Wang W, Xu Y, Gao R, Lu R, Han K, Wu G, et al. Detection of SARS-CoV-2 in Different Types of Clinical Specimens. *Jama* 2020;323:1843-4.
8. Zou L, Ruan F, Huang M, Liang L, Huang H, Hong Z, et al. SARS-CoV-2 Viral Load in Upper Respiratory Specimens of Infected Patients. *N Engl J Med* 2020;382:1177-9.
9. European Union Law. Directive 98/79/EC of the European Parliament and of the council of 27 October 1998 on in vitro medical devices. 1998. p. p.1-37.
10. Verheyen J, Kaiser R, Bozic M, Timmen-Wego M, Maier BK, Kessler HH. Extraction of viral nucleic acids: comparison of five automated nucleic acid extraction platforms. *J Clin Virol* 2012;54:255-9.
11. Torres I, Albert E, Navarro D. Pooling of nasopharyngeal swab specimens for SARS-CoV-2 detection by RT-PCR. *J Med Virol* 2020.
12. European Commission. Guidelines on COVID-19 in vitro diagnostic tests and their performance. Brussel: Communication from the commission. 2020.
<https://ec.europa.eu/info/sites/info/files/testingkitscommunication.pdf>

Conflict of interest

None.

Funding

None.

