

Article

Analysis of Glabrous Canary Seeds by ELISA, Mass Spectrometry and Western Blotting for the Absence of Cross-reactivity with Major Plant Food Allergens

Joyce Irene Boye, Allaoua Achouri, Nancy Raymond, Chantal Cleroux,
Dorcas Weber, Terry Koerner, Pierre Hucl, and Carol Ann Patterson

J. Agric. Food Chem., **Just Accepted Manuscript** • DOI: 10.1021/jf305500t • Publication Date (Web): 27 May 2013

Downloaded from <http://pubs.acs.org> on June 4, 2013

Just Accepted

"Just Accepted" manuscripts have been peer-reviewed and accepted for publication. They are posted online prior to technical editing, formatting for publication and author proofing. The American Chemical Society provides "Just Accepted" as a free service to the research community to expedite the dissemination of scientific material as soon as possible after acceptance. "Just Accepted" manuscripts appear in full in PDF format accompanied by an HTML abstract. "Just Accepted" manuscripts have been fully peer reviewed, but should not be considered the official version of record. They are accessible to all readers and citable by the Digital Object Identifier (DOI®). "Just Accepted" is an optional service offered to authors. Therefore, the "Just Accepted" Web site may not include all articles that will be published in the journal. After a manuscript is technically edited and formatted, it will be removed from the "Just Accepted" Web site and published as an ASAP article. Note that technical editing may introduce minor changes to the manuscript text and/or graphics which could affect content, and all legal disclaimers and ethical guidelines that apply to the journal pertain. ACS cannot be held responsible for errors or consequences arising from the use of information contained in these "Just Accepted" manuscripts.



ACS Publications
High quality. High impact.

Journal of Agricultural and Food Chemistry is published by the American Chemical Society, 1155 Sixteenth Street N.W., Washington, DC 20036
Published by American Chemical Society. Copyright © American Chemical Society. However, no copyright claim is made to original U.S. Government works, or works produced by employees of any Commonwealth realm Crown government in the course of their duties.

Analysis of Glabrous Canary Seeds by ELISA, Mass Spectrometry and Western Blotting for the Absence of Cross-reactivity with Major Plant Food Allergens

Joyce Irene Boye^{1*}, Allaoua Achouri¹, Nancy Raymond¹, Chantal Cleroux², Dorcas Weber², Terence B. Koerner², Pierre Hucl³, Carol Ann Patterson⁴

¹Food Research and Development Centre, Agriculture and Agri-Food Canada, 3600 Casavant Boulevard West, Saint-Hyacinthe, Qc, CANADA, ²Food Research Division, Bureau of Chemical Safety, Health Canada, Ottawa, CANADA, ³Crop Development Centre, University of Saskatchewan, Saskatoon, Saskatchewan, S7N 5A8, ⁴The Pathfinders Research and Management Limited, Saskatoon, Saskatchewan S7N0S5, Canada.

***Corresponding author:**

Tel: +1 450 768 3232

Fax: +1 450 773 8461

E-mail address: joyce.boyce@agr.gc.ca

ABSTRACT: Glabrous (hairless) canary seed belong to the *Poaceae* (*Gramineae*) family and could serve as an alternative source of gluten-free cereal grain. In this study, allergenic cross-reactivities between hairless, dehulled canary seeds (*Phalaris canariensis*) and major allergenic proteins from gluten, soy, peanuts, tree nuts, sesame and mustard were studied using commercial enzyme-linked immune sorbent assay (ELISA) kits specific for these target allergens. Mass spectrometry (MS) and immunoblotting were further used to assess for the presence of gluten specific protein fragments. MS results revealed the likely presence of proteins homologous with rice, oat, corn, carrot, tomato, radish, beet, and chickpea. However, no presence of celiac related gluten fragments from wheat, rye, barley or their derivatives were found. Immunoblotting studies yielded negative results further confirming the absence of gluten in the canary seed samples tested. No cross-reactivities were detected between canary seeds and almond, hazelnut, mustard, peanut, sesame, soy, walnut and gluten using ELISA.

KEYWORDS: Glabrous canary seed, gluten-free, gluten sensitive enteropathy, celiac disease

35 INTRODUCTION

36 Gluten sensitive (GS) enteropathy or celiac disease (CD) is an immune mediated
37 food sensitivity which affects about 0.3 - 1% of the population.^{1,2} CD is considered an
38 autoimmune disease because of the presence of tTG-targeting autoantibodies in the
39 serum and the intestinal mucosa of patients.³ Symptoms of CD include aneamia,
40 gastrointestinal disorders, stunted growth, dermatitis, endocrine disorders, amongst
41 many others.⁴ In practical terms the disease can be considered a spectrum disorder with
42 some individuals showing no symptoms whereas others experience extreme
43 debilitating symptoms. The offending proteins in CD is found in the gliadin fraction of
44 gluten (i.e., prolamin in wheat, secalin in rye and hordein in barley).⁵⁻⁷

45 To avoid symptoms and disease exacerbation, GS individuals must avoid all
46 sources of gluten in the diet. Gluten avoidance can however be a major challenge as
47 cereals containing gluten (e.g., wheat, barley, rye, spelt, kamut, triticale and their
48 derivatives) are frequently used in food preparation (e.g., as flours in bakery products,
49 sauces and seasonings, anti-caking agents, spices, icing sugar, etc). This limits choice
50 and further puts the nutritional status of CD patients at risk.

51 Growing awareness of CD and its impacts on health has garnered interest to
52 identify non-gluten containing cereal that could be safe for GS individuals. Non-gluten
53 containing cereals frequently used in product formulation include rice, teff, quinoa,

millet, buckwheat, sorghum, maize and amaranth. As the environmental conditions for growing these grains are variable, availability of regular supplies is not always assured.

An alternative grain that may potentially be considered for GS individuals is glabrous canary seed (*P. canariensis* L.). Annual canary seed belongs to the *Poaceae* (*Gramineae*) family, the Pooideae subfamily and the Agrostideae tribe. It is different from reed canarygrass (*Phalaris arundinacea* L.) which is a perennial forage grass of the same family and genus.⁷ There are over 15 species in the *Phalaris* genus, and of the annual species of the genus, *Phalaris canariensis* is the only one grown as a grain crop.⁸ The glabrous (hairless) trait and yellow coloured grain trait were derived by mutagenesis of the certified pubescent (hairy) parent *P. canariensis*, cultivar Keet.⁹

The chemical composition of canary seed shows that it has potential as a food crop. Glabrous canary seed contains an average of 24% protein, 8% crude fat, 56% starch and 7% total dietary fibre.¹⁰ Phytate levels are about twice that of wheat, and levels of trypsin inhibitor and alpha-amylase inhibitor are similar to those found in other cereal grains.¹¹ There are also smaller amounts of soluble sugars and ash in this grain.¹² Indeed, canary seed is a richer source of most required minerals relative to wheat.¹⁰ Baking tests show that bread made with 100% hairless canary seed flour is significantly lower in loaf volume and crust and crumb colour than wheat bread. Bread with loaf volume, specific volume and crust colour comparable to those of the wheat control was achieved by using up to 25% of hairless canary seed or 15% of roasted

canary seed flour, thus demonstrating its potential to replace a percentage of wheat flour in products or 100% in products such as cookies, cakes and flat breads for GS individuals.

Other more recent studies also show that hairless canary seed has significantly higher levels of phosphorus, sulphur, magnesium, calcium, iron, manganese and zinc than wheat and oat. Thus, from a nutritional perspective, glabrous canary seed could be a good source of macronutrients such as protein and fiber as well as a source of some important micronutrients including phenolic compounds and carotenoids.¹³⁻¹⁵

Reports of the pollen of perennial forage canarygrass (e.g., *Phalaris aquatica*, *Phalaris arundinacea*) as a major environmental allergen and incidents of allergic reactions to canaryseed on inhalation during handling have been cited in the literature. Using IgE antibodies from sera of 24 grass-pollen-allergic subjects, seventeen allergenic fractions of canarygrass (*Phalaris aquatica*) pollen were identified ranging in molecular mass from 14 to 100 kDa. A 34-kDa protein fraction was found to have the highest frequency of IgE binding (77%) and was tentatively designated as Pha a I.¹⁶ Another study reported immediate rhinoconjunctivitis, asthma, and contact urticaria while handling bird food caused by black seed (*G. abyssinica*) used as ingredient in canary food.¹⁷

In order to establish the safety of canary seeds for human consumption from a food allergy perspective and further allow the mitigation of potential risks which could arise with their consumption, studies are required to confirm the absence of cross-reactivity with major allergenic proteins. The specific objective of this study, therefore, was to assess the safety of glabrous canary seed from an allergy perspective by analysing for the presence of allergenic cross-reactivities between glabrous canary seed and allergenic proteins in the following foods: gluten, soy, peanuts, tree nuts, sesame and mustard. Mass spectrometry (MS), bioinformatics and immunoblotting were used to assess the homology of glabrous canary seed proteins with those from known gluten sources (wheat, barley and rye). Additionally, commercially available ELISA kits were used for the detection of cross-reactivities with each target food allergen as indicated in Tables 1 and 2.

MATERIALS AND METHODS

Materials and Chemicals. Hairless canary seed (*P. canariensis* L.) cultivars, CDC Maria (brown coloured grain), and two yellow coloured grain cultivars, C05041 and C05091, were used in this study. They were grown in three-replicate randomized complete block experiments in three locations in Saskatchewan, Canada over a two year period (2007 & 2008). Harvested samples from each of the test plots were dehulled and hand cleaned to remove any broken seeds or foreign material. Each sample was a blend

of seeds from the same variety grown on different plots. Further cleaning was done by spraying the surface of the seeds with 70% ethanol, gentle cleaning using paper towels, followed by air-drying at room temperature in a fume hood. The samples were ground in a WARING commercial blender 7011 (Model 31BL92, Waring products division, New Hartford, Connecticut, 06057, U.S.A) and stored at -20°C until analyzed. Oat, millet, quinoa, teff, sorghum and buckwheat (marketed as gluten-free) and Canadian western red spring wheat were also obtained through local suppliers, cleaned and milled with a coffee grinder and similarly stored for western blotting analysis.

Serum samples of 10 patients with wheat allergy were purchased from PlasmaLab International (Everett, WA, U.S.A.). The sera were portioned in aliquots and stored at -80°C until use. Mouse Anti-Human IgE (monoclonal) was purchased from SouthernBiotech (Birmingham, AL, USA). Wheat gluten and β -lactoglobulin polyclonal antibodies and goat anti-rabbit secondary antibodies conjugated to horse radish peroxidase (HRP) were purchased from Cederlane (Burlington, Ontario, Canada) and Abcam Inc. (Cambridge, MA, USA), respectively. Rabbit polyclonal anti-sesame antibodies were donated by Health Canada Proteomic and Allergen Research Unit.

Dithiothreitol (DTT), iodoacetamide, [Glu1]-Fibrinopeptide B, ammonium bicarbonate (NH_4HCO_3), acetonitrile (ACN) and formic acid (FA) were purchased from Sigma-Aldrich (Oakville, ON) and used without purification. Promega sequencing grade trypsin was purchased from Fisher Scientific (Nepean, ON) and ethanol was

purchased from Commercial Alcohols (Brampton ON). Water was purified with resistivity $\geq 18 \text{ M}\Omega$ using a Millipore Milli-Q Synthesis system from Fisher Scientific (Nepean, ON). The extraction solution used in these experiments consisted of 60% ethanol and buffer (NH_4HCO_3 , 100 mM, pH 8.5 with 1M Urea). Molecular mass calibration kit (10 to 250 kDa) from Bio-Rad (Bio-Rad Laboratories, Hercules, CA, USA) was used as molecular markers.

ELISA kits for selected priority allergens including standards were purchased from commercial suppliers as listed in Table 1. Information for each kit and the websites of each kit manufacturer for full protocols are provided in Tables 1 and 2. All other reagents and chemicals used were of analytical grade.

Sample Preparation for Mass Spectrometric Analysis. A sample of finely ground canary seed powder (5 g) was dispersed in the extraction solution (5 mL) by shaking on a vortex for 30 seconds. An aliquot of the supernatant (1mL) was diluted and the protein content determined using a BCA test from Fisher Scientific (Nepean, ON). Further dilutions were made to produce a standard concentration ($25 \mu\text{g}/100 \mu\text{L}$) for the digestion and subsequent mass spectrometry analysis. Dithiothreitol (DTT; 10 mM in 100 mM ammonium bicarbonate) was added for reduction, while iodoacetamide (IA; 55 mM in 100 mM ammonium bicarbonate) was added for alkylation. Samples were then digested with trypsin, and finally, 1% formic acid (FA) was added to the mixture prior to MS analysis.

153 **LC-MS/MS.** The LC-MS/MS system used was a Waters CapLC system coupled
154 to a Waters quadrupole-time of flight mass spectrometer (QTOF2) with nanolockspray
155 (Waters, Milford, MA). The mass spectrometer was operated in nanoelectrospray
156 positive ion detection mode with a spray voltage of 3.5 kV at a resolution of 10000
157 (FWHH). The LC system was configured with a binary solvent manager to deliver the
158 gradient system and an auxiliary solvent manager to deliver the reference mass
159 standard solution for accurate mass corrections. The reference mass standard consisted
160 of a solution of [Glu1]-Fibrinopeptide B in 50% aqueous methanol, which was infused
161 into the source at a flow rate of 0.200 $\mu\text{L}/\text{min}$ and sampled every 10 seconds for one
162 second. The separation system consisted of a trap column (Symmetry C18 5 μm 180 μm
163 x 20 mm, Waters) and an analytical column (BEH 130 C18 1.7 μm 100 μm x 100 mm,
164 Waters). Solvent A consisted of H_2O with 0.1% FA and solvent B consisted of ACN
165 with 0.1% FA. In a typical experiment digested canary seed protein samples were
166 concentrated and desalted by injection (1 μL) onto the trap column for 2 min at a flow
167 rate of 15.0 $\mu\text{L}/\text{min}$ using 99% solvent A. The samples were then diverted to the
168 analytical column and eluted at 400 nL/min. The elution program started with 99% A
169 for 1 min followed by a gradient from 1% to 50% B in 50 min. The mass spectrometer
170 was operated in positive ion mode (resolution >10000 FWHH) and the acquired data
171 collected using MassLynx v4.1. Data was collected in both MS survey scan mode and an
172 automated data directed analysis (DDA) mode. MS survey scan data was acquired in

continuum mode from m/z 400 to 1600 with a scan time of 0.60 seconds with an interscan delay of 0.02 seconds.

Data collected using data-directed analysis (DDA) software was submitted directly to an in-house MASCOT server with software version 2.3.01. The tandem MS information was searched against the National Center for Biotechnology Information (NCBI) database (version 20100115) with up to 2 missed cleavages, a fixed carbamidomethylation modification at cysteine and oxidation of methionine, and deamidation of asparagine and glutamine were specified as variable modifications. The peptide mass tolerance and the fragment mass tolerance were both set to 0.1Da. The mass and the retention time data are collected for peptides that produced a significant identification or extensive homology (probability-based MOWSE score ≥ 38).

Identification of Protein Bands as Allergens by LC-MS/MS. The protein in-gel digestion and mass spectrometry experiments were performed by the Proteomics platform of the Eastern Quebec Genomics Center, Quebec, Canada. Detailed experimental parameters for the tryptic digestion, the mass spectrometry conditions and data analysis were previously reported.¹⁸ Scaffold (Scaffold_3_00_07, Proteome Software Inc., Portland, OR) was used to validate MS/MS based peptide and protein identifications. Peptide identifications were accepted if they could be established at greater than 95.0% probability as specified by the Peptide Prophet algorithm.¹⁹ Protein identifications were accepted if they could be established at greater than 95.0%

193 probability and contained at least 2 identified peptides. Protein probabilities were
194 assigned by the Protein Prophet algorithm.²⁰

195 **SDS-PAGE and Western Blot Analysis.** SDS-PAGE and immuno-blotting were
196 used for the determination of allergenic proteins. SDS-PAGE was performed according
197 to Laemmli²¹ using 10-20% precast gel (Bio-Rad, St. Louis, Mo., U.S.A.). Protein samples
198 (2.5mg/mL) were dispersed in Laemmli sample buffer (Bio-Rad) plus 5% β -
199 mercaptoethanol and vortexed for 15 min. The dispersions were heated at 100 °C for 5
200 min, and centrifuged before loading. 10 μ L of the supernatant of each sample was
201 loaded in a well and the gel was run at constant voltage of 150V for approximately 90
202 min at room temperature. After electrophoresis, the separated proteins were stained
203 with Coomassie Brilliant Blue R-250 and scanned.

204 For immuno-blotting, the separated proteins were subsequently transferred to a
205 polyvinylidene difluoride (PVDF) membrane (Millipore, 0.45 mm) using a Criterion
206 Blotter cell (Bio-Rad). The transfer was done for 1 h with a constant current of 100 V in
207 25 mM Tris-HCl (pH 8.3), 192 mM glycine, 20% methanol, following the method of
208 Towbin et al.²² Blots were blocked in phosphate-buffered saline plus 0.1% Tween 20
209 containing 5% non-fat milk (NFM) powder or for 1 h at room temperature, and then
210 the PVDF was washed three times with 0.1% Tween 20 in TBS and incubated overnight
211 at 4 °C with a pooled sera from wheat-allergic individuals (diluted 1:30 with 0.1%
212 Tween 20 in TBS containing 1% BSA protease free) or with gluten-IgG (dilution 1:1000),

213 β -lactoglobulin-IgG (dilution 1:10000) or sesame-IgG (dilution 1:30000) The membrane
214 was washed and incubated for an additional 1 h with a goat anti-human IgE polyclonal
215 (HRP) conjugate used as second antibody (diluted 1:5000) (KPL, Inc. Gaithersburg,
216 Maryland, U.S.A.). Incubations were performed with gentle rocking. After washing, the
217 blot was incubated in the Bio-Rad amplification reagent (BAR), washed, incubated in
218 streptavidin-HRP, and washed again before incubation with the Opti-4CN substrate kit
219 (Bio-Rad) to visualize the immunoreactive bands. To improve binding sensitivity and
220 reduce non-specific binding, a mouse anti-human IgE monoclonal antibody (dilution
221 1:2000) as well as diluting all antibodies in the blocker buffer (i.e. 3% non-fat dry milk)
222 were also tested in comparison with the originally used goat anti-human IgE and 1%
223 BSA diluting buffer.

224 Both the Coomassie-stained gels and immunostained membranes were scanned
225 using an Epson Image Scanner III (Seiko Epson, Corporation, Nagano, Japan), and the
226 generated files were analyzed with ImageQuant TL software (Amersham Biosciences)
227 using the Kaleidoscope Precision Plus Protein Standards from Bio-Rad.

228 For the immunoblotting studies, the presence of a band in the immunoblot
229 indicated the likely presence of an IgE antibody-antigen binding interaction.

230 **Enzyme Linked Immuno Sorbent Assay (ELISA).** ELISA analysis was carried
231 out as per the instructions provided by the kit manufacturers (Tables 1 and 2). In
232 general, glabrous canary seed were ground into flour and extracted as per the kit

instructions. Samples were dispersed in the extraction solutions provided with the ELISA kit and were continuously stirred while maintained at the temperature specified in the kit instructions and for the periods specified. The extracts were then filtered and tested according to the kit instructions. The protocol used varied slightly for the different test kit depending on whether the ELISA assay was direct, sandwich or competitive. When cross-contamination was suspected, samples were cleaned with 70% ethanol and the extractions were repeated. To augment sensitivity for further confirmation, the proposed amounts indicated on the kit instruction were tripled and the extractions were repeated for some samples. Optical densities (OD) were measured at 450-650 nm depending on the kit instructions using a Bench Mark Plus microplate reader (Bio-Rad Laboratories Ltd., Mississauga, Ont., Canada) controlled by a Microplate Manager Software (Bio-Rad Laboratories Ltd., Mississauga, Ont., Canada). Due to reported variability in ELISA results from different test kits, two to three commercial test kits from different companies were used for each targeted allergen and extractions were done in triplicate for each kits and each extract was analysed in triplicate.

Quantitative determination of allergen content for ELISA was made using the standards provided in each kit and, when available, using the accompanying software provided by the kit manufacturer, otherwise the Excel software was used to calculate the linear/exponential regression slopes of the standards and to determine the allergen

content. Values below the limit of detection (LOD) or limit of quantitation (LOQ) were reported as such indicating the absence of the target allergen. Values above the LOD or LOQ were reported as averages of all determinations. This indicated the likely presence of a target allergen. Data were statistically evaluated by 1-way analysis of variance (ANOVA) using the PRISM software, version 3.02 (Graph Pad Software, Inc. San Diego, Calif., U.S.A.). Significant differences between means were determined by the Tukey's Multiple Comparison Test procedure at the 5% significance level.

RESULTS AND DISCUSSION

Mass Spectrometry. The mass spectrometry MASCOT database search scores for protein hits were cut off at 38 which represent significant homology between the test sample and identified protein matches. Tables 3 presents a list of the proteins identified from the database in the three canary seed samples which were mostly homologous with rice, oats, corn, carrot, tomato, radish, beet, and chickpea proteins. The search showed no significant homology to celiac related gluten fragments from wheat, rye, barley or their derivatives in all three canary seed samples. For CDC Maria (Table 3), three hits were obtained indicating the likely presence of protein disulfide-isomerase (wheat), Em protein H5 (wheat) and cytosolic glyceraldehyde-3-phosphate dehydrogenase (barley) or proteins having similar homology. One hit suggesting the likely presence of cytosolic glyceraldehyde-3-phosphate dehydrogenase (barley) or a

273 similar protein was found for C05041 (Table 3). Protein disulfide-isomerase, with a
274 molecular mass of 56,533 Da, is an enzyme in the endoplasmic reticulum in eukaryotes
275 that catalyzes the formation and breakage of disulfide bonds between cysteine residues
276 within proteins as they fold.²³ Em protein H5 (molecular mass, 10,060 Da) is a member
277 of the small hydrophilic plant seed protein family. Cytosolic glyceraldehyde-3-
278 phosphate dehydrogenase (molecular mass, 33,236 Da) belongs to the glyceraldehyde-3-
279 phosphate dehydrogenase family.

280 Gluten epitopes provoking celiac disease typically originate from the gliadin and
281 glutenin fractions and contain high amounts of glutamine and proline amino acid
282 residues and the signature “QP” amino acid sequence.²⁴ Examples of key epitopes
283 reported in the literature include IIQPEQPAQ, FPQQPQQPYYPQQP,
284 FSQPQQQFPQPQ, LQPQQPFPQQPQQPYYPQQPQ, PQQPFPQPQQQFPQPQQPQQ.²⁵
285 The amino acid sequences of the three protein hits (i.e., protein disulfide-isomerase
286 (wheat), Em protein H5 (wheat) and cytosolic glyceraldehyde-3-phosphate
287 dehydrogenase (barley) (Figure 1C), do not show any “QP” sequence suggesting little
288 likelihood of them containing a celiac provoking epitope. Overall, the mass
289 spectrometry results suggest either cross contact or homology between canary seed
290 proteins and some proteins from rice, oats, corn, carrot, tomato, radish, beet, and
291 chickpea proteins but not with gluten.

Immunoblotting. The SDS-PAGE protein profiles of the three canary seed samples are shown in Figure 1A. All three seeds showed similar electrophoretic profiles with the presence of protein bands ranging in molecular mass from ~ 10,000 to 100,000 Da. The most prominent band had a molecular mass of ~ 20,000 – 25,000 Da.

Reactivities of proteins separated by SDS-PAGE were analyzed using pooled sera from 10 wheat allergenic individuals (Figure 1B, C, and D). The amplified Opti-4CN detection kit used to perform colorimetric blotting requires the use of 1% BSA as diluent buffer for primary and secondary antibodies as well as the substrate reagents. The preliminary immunoblot of the three canary seed protein extracts using bovine serum albumin as antibody diluent buffer revealed strong binding with many of the canary seed protein bands (Figure 1B). These results were intriguing as the mass spectrometry results showed the absence of celiac provoking gluten proteins. Binding, however, could sometimes result from non-specific binding involving one of the detecting reagents used during the assay (e.g., BSA). Indeed, blotting of canary seed with only the secondary antibodies (no primary human sera used) resulted in the binding of some of the canary seed proteins, confirming the presence of non-specific binding (data not shown).

Thus, the blotting was repeated this time using non-fat dry milk (NFM), Pierce® protein-free blocking and superbloc buffers as buffer diluent and blocking agents, respectively. Negative controls were also tested to reveal any non-specific binding in

the absence of human sera. The NFM blocking and diluent buffer showed good binding to all potential sites of nonspecific interaction, eliminating background altogether without altering or obscuring the epitope for antibody binding (Figure 1D). Furthermore, as previously reported by other workers,²⁶ the use of a monoclonal secondary antibody (mouse anti-human IgE) gave better results with higher sensitivity than the polyclonal secondary antibody (goat anti-human IgE) (Figures 1B and 1C).

Optimal conditions for the immunoblot analysis of the canary seed with wheat human sera were thus obtained with the monoclonal secondary antibody and NFM as blocker (Figure 1D). As shown, the results revealed that some canary seed proteins were still recognised by the pooled sera of wheat allergic patients. Two strong antibody binding proteins with molecular mass (MM) ranging between 50 to 75 kDa were observed in addition to other polypeptides having MMs of 37-50 kDa and 15 kDa. The three canary seed composites showed similar antibody-binding patterns.

Aalberse et al.²⁷ reported homologous allergens from phylogenetically related grasses. To determine if binding will be observed with gluten specific antibodies, the work was repeated with polyclonal rabbit gluten IgG antibodies raised specifically against wheat gluten protein (immunogen). In addition, pooled sera of 7 individuals allergic to sesame as well as anti β -lactoglobulin antibody were tested as negative controls. As shown in Figure 2, despite the wheat sample used as positive control, no

bands were observed in all three immune blots suggesting the absence of gluten specific proteins in the three canary seed samples.

Immunoblotting of Other Gluten-Free Cereals. To verify if binding will occur with other gluten-free cereals, the SDS-PAGE and blotting was performed on oat, millet, quinoa, teff, sorghum and buckwheat as shown in Figure 3. Canadian western red spring wheat was used as positive control. The SDS-PAGE results (Figure 3A) showed major differences in the electrophoretic profiles of the cereals. This was expected as the cereals belong to different plant families. As observed for canary seeds, when bovine serum albumin was used as blocking agent the pooled wheat sera recognized practically all the different polypeptide bands that were clearly visible in the SDS-PAGE profile of the other gluten-free cereals as well as some that were not previously evident (Figure 3B). Blocking with non-fat dry milk instead of bovine serum albumin revealed a different pattern with only few bands recognised (Figure 3C). Prominent bands at approximately 50 to 75 kDa showed similar strong antibody-binding. Faint bands (low binding) were also observed for some proteins with MM above 37 kDa and at ~15 kDa. The western blotting was repeated using rabbit polyclonal gluten antibodies with non-fat dry milk as blocking agent. The immunoblot (Figure 3D) revealed strong binding to many of the wheat proteins and some proteins in oat, millet, quinoa, teff, and to a lower extent with sorghum and buckwheat proteins, which could be either due to cross-reactivities or cross-contamination of the grains with gluten proteins.

In previous studies using radioallergosorbent testing (RAST) of sera from subjects sensitized to wheat and rye flour, a significant reaction was found with seed extracts of 12 cereals including wheat, durum wheat, triticale, cereal rye, barley, rye grass, oats, canarygrass, rice, maize, sorghum and Johnson grass.²⁸

Identification of Major Allergens (excised bands) by Mass Spectrometry.

To confirm the identity of the predominant protein components recognised by antibodies in the wheat sera, electrophoresis of wheat, gluten-free cereals and canary seeds was run again and the bands showing antibody-antigen binding during immunoblotting (Figure 4) were excised and further analyzed by LC/ESI-MS/MS. The list of proteins identified is presented in Table 4. Specifically, for the canary seeds, the allergenic protein bands (as shown in Figure 4B) were not identified by the MS technique as belonging to the *Phalaris canariensis* family. To date, very few proteins from canaryseeds have been sequenced. However, using Mascot search, the tryptic peptides from the canary bands (C1-C7) were compared with peptides of known proteins listed in the National Center for Biotechnology Information (NCBI), and the results revealed significant sequence homology to several cereals in particular with rice (20% coverage), oat (16%), barley (14%), sorghum (10%), maize (5%) and granule-bound starch synthase 1 (GBSS1) from wheat (4%).

ELISA. ELISA results of the canary seeds for the different allergen kits tested and using the sample size indicated in the kit instructions are provided in Table 5. All the

results were below the LOD/LOQ, except for two test kits (i.e., ELISA System kit for mustard and Neogen Veratox for gluten and mustard). The ELISA System mustard test kit gave significantly ($p < 0.05$) higher values for mustard ranging between 0.7 – 13 ppm for the different canary seed composites. The LOQ for this kit was 1.0 ppm. The Neogen Veratox test kits also gave value of 5.4 ppm (close to LOQ) for gluten and 36.5 ppm for mustard for CDC Maria which was above the LOQ. A significant difference ($p < 0.05$) was observed between the mustard levels determined by the ELISA System kit and Neogen Veratox particularly for the CDC Maria variety. For the two other composites (C05041 and C05091) the Neogen Veratox results were below the LOQ of the kits (5 ppm for gluten and 2.5 ppm for mustard).

As the results for mustard and gluten using the sample size indicated in the kit instructions were above the LOQ (Table 5), the concentration requested by the kit instruction was tripled and the analysis was repeated. This was done intentionally in order to augment the sensitivity of the detection. Additionally, the canary seed groats were cleaned by washing with 70% ethanol, air-dried and reanalysed as per the kit instructions using one and three times the sample size recommended for extraction. The ELISA results obtained (Tables 6 and 7) were below the LOD/LOQ for all the kits re-tested indicating the absence of cross-reactivity or cross-contamination with the targeted allergens. The results further suggest that the earlier positive detection was

390 likely due to localised cross-contamination with some gluten and mustard containing
391 material in the sample or a false positive.

392 A new Veratox kit for gluten which uses the Codex Alimentarius recommended
393 “R5” gliadin antibody was used to test for the presence of gluten in the three different
394 canary seed composites, and the results as shown in Table 7, confirmed the absence of
395 cross-reactivity or gluten cross-contamination with the analyzed canary seed
396 composites. Samples of the same groats were sent to the company for testing and their
397 results (not shown) further confirmed the absence of gluten.

398 In summary, glabrous canary seed could be a good alternative gluten-free cereal.
399 CODEX²⁹ guidelines recommend that foods containing less than 20 ppm gluten may be
400 labeled as gluten free. All three techniques used in this study were negative for gluten,
401 thus glabrous canary seeds could be safely labelled as gluten free. As glabrous canary
402 seed is likely to be grown in close proximity to other gluten containing crops such as
403 wheat, rye and barley, the risk of cross-contamination during production, harvesting
404 and processing can be high if appropriate allergen management practices are not used.
405 Food processors interested in using this grain in gluten-free foods would need to ensure
406 that the grains are appropriately managed from production, transportation through
407 primary processing in order to ensure that they remain “gluten-free”. Interestingly, the
408 immunoblotting studies showed that non-gluten containing cereals such as buckwheat,
409 sorghum, teff, millet, quinoa and oat contain proteins that are recognised by antibodies

in wheat allergic sera and rabbit anti-gluten IgG raised against wheat gluten proteins. As IgE mediated wheat allergy is different from GSE, the mechanisms underlying these binding interactions require further studies.

ACKNOWLEDGEMENTS

This project was funded by the CanaryseedDevelopment Commission of Saskatchewan. In-kind support from Agriculture & Agri-Food Canada is gratefully acknowledged. We would also like to thank Dr. Sahul H. Rajamohamed for his assistance with the preparation of the references for this report.

REFERENCES

- (1) Chin, L. C.; Jones, M. K.; Kingham, J. G. C. Celiac disease and autoimmune thyroid disease. *Clin. Med. Res.* **2007**, 5, 184-192.
- (2) Di Sabatino, A.; Corazza, G. R. Celiac disease. *The Lancet* **2009**, 373, 1480-1493.
- (3) Caputo, I.; Lepretti, M.; Martucciello, S.; Esposito, C. Enzymatic Strategies to Detoxify Gluten: Implications for Celiac Disease. *Enz. Research.* **2010**, 5, 1-9.
- (4) Kagnoff, M. F. Celiac disease: Pathogenesis of a model immunogenetic disease. *J. Clin. Invest.* **2007**, 117, 41-49.

- 428 (5) Molberg, Ø.; Flæte, N. S.; Jensen, T.; Lundin, K. E. A.; Arentz-Hansen, H.;
429 Anderson, O. D.; Kjersti Uhlen, A.; Sollid, L. M. Intestinal T-cell responses to high-
430 molecular-weight glutenins in celiac disease. *Gastroenterology* **2003**, 125, 337-344.
- 431 (6) Dewar, D. H.; Amato, M.; Ellis, H. J.; Pollock, E. L.; Gonzalez-Cinca, N.; Wieser,
432 H.; Ciclitira, P. J. The toxicity of high molecular weight glutenin subunits of wheat to
433 patients with coeliac disease. *Eur. J. Gastroenterol. Hepatol.* **2006**, 18, 483-491.
- 434 (7) Wieser, H. Chemistry of gluten proteins. *Food Microbiol.* **2007**, 24, 115-119.
- 435 (8) Putnam, D. H.; Miller, P. R.; Hucl, P. Potential for production and utilization of
436 annual canarygrass. *Cereal Foods World* **1996**, 41, 75-83.
- 437 (9) Hucl, P.; Matus-Cadiz, M.; Vandenberg, A.; Sosulski, F. W.; Abdel-Aal, E. S. M.;
438 Hughes, G. R.; Slinkard, A. E. CDC Maria annual canarygrass. *Can. J. Plant Sci.* **2001**, 81,
439 115-116.
- 440 (10) Abdel-Aal, E. S. M.; Hucl, P.; Miller, S. S.; Patterson, C. A. and Gray, D.
441 Microstructure and nutrient composition of hairless canary seed and its potential as a
442 blending flour for food use. *Food Chem.* **2011a**, 125, 410-416.
- 443 (11) Abdel-Aal, E. S. M.; Hucl, P.; Patterson, C. A.; Grey, D. Phytochemicals and
444 heavy metals content of hairless canary seed: a variety developed for food use. *LWT-*
445 *Food Sci. Technol.* **2011b**, 44, 904-910.

- 446 (12) Abdel-Aal, E. S. M.; Hucl, P.; Sosulski, F. W. Structural and compositional
447 characteristics of canary seed (*Phalaris canariensis* L.). *J. Agric. Food Chem.* **1997b**, *45*,
448 3049-3055.
- 449 (13) Abdel-Aal, E.S.M., Hucl, P. Hairless canary seed: a potential food crop, In *Specialty*
450 *grains for food and feed*, Ed. Abdel-Aal, E., Wood, P., American Association of Cereal
451 Chemists: St. Paul, Minnesota, **2005**. 203-221.
- 452 (14) Li, W.; Qiu, Y.; Patterson, C. A.; Beta, T. The analysis of phenolic constituents in
453 glabrous canary seed groats. *Food Chem.* **2011**, *127*, 10-20.
- 454 (15) Li, W.; Beta, T. An evaluation of carotenoid levels and composition of glabrous
455 canary seed. *Food Chem.* **2012**, *133*, 782-786.
- 456 (16) Suphioglu, C.; Singh, M. B.; Simpson, R. J.; Ward, L. D.; Knox, R. B. Identification
457 of canary grass (*Phalaris aquatica*) pollen allergens by immunoblotting: IgE and IgG
458 antibody-binding studies. *Allergy* **1993**, *48*, 273-281.
- 459 (17) Rodríguez, B.; Rodríguez, A.; de Barrio, M.; Tornero, P.; Baeza, M. L. Asthma
460 induced by canary food mix. *Allergy Asthma Proc.* **2003**, *24*, 265-268.
- 461 (18) Rioux, G.; Majeau, N.; Leclerc, D. Mapping the surface-exposed regions of
462 papaya mosaic virus nanoparticles. *FEBS Journal* **2012**, *279*, 2004-2011.
- 463 (19) Keller, A.; Nesvizhskii, A. I.; Kolker, E.; Aebersold, R. Empirical statistical model
464 to estimate the accuracy of peptide identifications made by MS/MS and database search.
465 *Anal. Chem.* **2002**, *74*, 5383-5392.

- 466 (20) Nesvizhskii, A. I.; Keller, A.; Kolker, E.; Aebersold, R. A statistical model for
467 identifying proteins by tandem mass spectrometry. *Anal. Chem.* **2003**, 75, 4646-4658.
- 468 (21) Laemmli, U. K. Cleavage of structural proteins during assembly of the head
469 bacteriophage T4. *Nature* **1970**, 227, 680-685.
- 470 (22) Towbin, H.; Gordon, J. Immunoblotting and dot blotting: current status and
471 outlook. *J. Immunol.* **1984**, 72, 313-340.
- 472 (23) Wilkinson, B.; Gilbert, H. F. Protein disulfide isomerase. *Biochim. Biophys. Acta*,
473 **2004**, 1699, 35-44.
- 474 (24) Osman, A. A.; Gué Nnel, T.; Dietl, A.; Uhlig, H. H.; Amin, M. Fleckenstein, B.;
475 Richter, T.; Mothes, T. B cell epitopes of gliadin. *Clin Exp Immunol.* **2000**, 121, 248-254.
- 476 (25) Qiao S. W.; Bergseng, E.; Molberg, Ø.; Jung, G.; Fleckenstein, B.; Sollid, L. M.
477 Refining the rules of gliadin T-cell epitope binding to the disease-associated DQ2
478 molecule in celiac disease: Importance of proline spacing and glutamine deamidation.
479 *J. Immunol.* **2005**, 175, 254-261
- 480 (26) Lipman, N. S.; Jackson, L. R.; Trudel, L. J.; Weis-Garcia, F. Monoclonal versus
481 polyclonal antibodies: Distinguishing characteristics, applications, and information
482 resources. *Inst. Lab. Anim. Res. J.* 2005, 46, 258-268.
- 483 (27) Aalberse, R. C.; Akkerdaas, J. H.; van Ree, R. Cross-reactivity of IgE antibodies to
484 allergens. *Allergy* **2001**, 56, 478-490.

485 (28) Baldo, B. A.; Krilis, S.; Wrigley, C. W. Hypersensitivity to inhaled flour allergens:
486 Comparison between cereals. *Allergy* **1980**, 35, 45-56.

487 (29) Codex Alimentarius Commission. Draft revised codex standard for foods for
488 special dietary use for persons intolerant to gluten. 2008 (ALINORM 08/31/26,
489 Appendix III). Codex Alimentarius Commission, 2008. Codex standard for foods for
490 special dietary use for persons intolerant to gluten. Codex Standard 118-1979.

491

492

493 LIST OF FIGURES

494 **Figure 1.** (A) SDS-PAGE patterns of canary seeds; (B) Immunoblot pattern of canary
495 seeds using pooled sera of 10 wheat allergic individuals. Antibodies and streptavidin
496 diluted in a non-fat dry milk buffer. Std - molecular weight standards; C1- CDC Maria;
497 C2- C05041; and C3- C05091.

498 **Figure 2.** Immunoblots of canary seeds (C1 – CDC Maria, C2 – C05041, C3 – C05091)
499 and W-wheat (positive control) after incubation with (A) Rabbit anti-gluten IgG
500 antibody, (B) Rabbit anti-beta-lactoglobulin IgG antibody, (C) Rabbit anti-sesame IgG
501 antibody. Antibodies and streptavidin diluted in a non-fat dry milk buffer. The
502 immunoblots were performed using Opti-4CN Colorimetric detection method.

503 **Figure 3.** (A) SDS-PAGE patterns of gluten and non-gluten containing cereals. (B, C)
504 Immunoblot patterns of gluten and non-gluten containing cereals. (B) Bound IgE
505 detected using pooled sera of 10 wheat allergic individuals. Antibodies and streptavidin
506 diluted in a non-fat dry milk buffer. (C) Immunoblot patterns bound using gluten
507 specific IgG antibodies (host, rabbit). Antibodies and streptavidin diluted in non-fat dry
508 milk. Std - molecular weight standards. The immunoblots were performed using Opti-
509 4CN Colorimetric detection method.

510 **Figure 4.** (A) Immunoblot patterns of gluten (wheat) and non-gluten containing cereals
511 (oat, millet, quinoa, teff, sorghum and buckwheat) and (B) Immunoblots of canary seeds

512 (C1 – CDC Maria, C2 – C05041, C3 – C05091). Bound IgE was detected using pooled
513 sera of 10 wheat allergic individuals.

514

515

516

517

518

519

520

521

522

Table 1. Enzyme Linked Immuno Sorbent Assay (ELISA) Commercial Test Kits Used for the Assessment of Allergen Cross Reactivity Between Glabrous Canary Seed and Selected Priority Allergens

Allergen	Company	Test kit	Result expression	Limit of detection/quantification (LOD)/(LOQ)	Range of quantification	Extraction method	Antibody
Almonds	Gen-Probe/Tepnel	Biokit	Almond	ppm 0.1 (LOD)	ppm 1 - 20	1 to 10 in Tris (hydroxymethyl) aminomethane / sodium chloride / gelatin	Polyclonal antibody specifically detects almond proteins
	Neogen	Veratox	Almond	2.5 (LOQ)	2.5 - 25	1 to 25 in PBS solution + additive	-
	R-Biopharm	Ridascreen	Almond	1.7 (LOD)	2.5 - 20	1 to 20 in buffer solution	Anti-almond antibodies
Gluten	Gen-Probe/Tepnel	Biokit	Gluten	1.0 (LOD)	3 - 50	1 to 10 in buffer solution with gluten extraction additive and 40% ethanol	Monoclonal (401/21) reactive to glutenin and gliadin
	Neogen	Veratox	Gliadin / Gluten	5.0 (LOQ)	5 - 50 gliadin	1:10 in 40% ethanol, then 1:40 in PBS + additive	-
	Neogen	Veratox	Gliadin R5	2.5 (LOQ)	2.5 - 40	Gliadin renaturing cocktail solution	-

	R-Biopharm	Ridascreen	Gliadin / Gluten	3.0 (gluten) (LOD)	5 ppb - 80 ppb gliadin	1 to 10 in 60% ethanol	Monoclonal R5 to gliadin-fractions
Hazelnut	Gen-Probe/Tepnel	Biokit	Hazelnut	0.1 (LOD)	1.6 - 32	1 to 10 in Tris (hydroxymethyl) aminomethane / sodium chloride / gelatin	Polyclonal antibody specifically detects hazelnut proteins
	Neogen	Veratox	Hazelnut	2.5 (LOQ)	2.5 - 25	1 to 25 in PBS solution + additive	-
	R-Biopharm	Ridascreen	Hazelnut	1.5 (LOD)	2.5 - 20	1 to 20 in buffer solution + 1g milk powder	-
Mustard	ELISA System	Mustard Seed Protein Residue	Mustard Seed Protein	1.0 (LOQ)	1 - 10	1 to 10 in extraction buffer solution	Anti-mustard seed protein antibodies
	Neogen	Veratox	Mustard	2.5 (LOQ)	2.5 - 25	1 to 25 in PBS solution + additive	-
	Sedium R&D	ELISA Kit	Mustard	0.06 (LOD)	0.5 - 15	1 to 10 in extraction buffer solution	-
Peanut	Gen-Probe/Tepnel	Biokit	Peanut	0.1 (LOD)	1 - 20	1 to 10 in Tris (hydroxymethyl) aminomethane / sodium chloride / gelatin	Polyclonal antibody specifically detects conarachin (Ara h 1)
	Neogen	Veratox	Peanut	2.5 (LOQ)	2.5 - 25	1 to 25 in PBS solution + additive	-

Sesame	R-Biopharm	Ridascreen	Peanut	1.5 (LOD)	2.5 - 20	1 to 20 in extraction buffer solution	Polyclonal antibody specifically detects peanut protein, including the peanut allergen Ara h 1
	ELISA System	Sesame Residue Detection	Sesame Seed	0.5 (LOQ)	0.5 - 5	1 to 10 in extraction buffer solution	Specific anti-sesame seed 2S-albumin antibodies
	Gen-Probe/Tepnel	Biokit	Sesame	0.1 (LOD)	6.125 - 100	1 to 10 in Tris (hydroxymethyl) aminomethane / sodium chloride / gelatin	Sesame antibodies
Soy	ELISA System	Soy Protein Residue Detection	Soy Flour Protein	2.5 (LOQ)	2.5 - 25	1 to 10 in extraction buffer solution	Specific anti-soy trypsin inhibitor and other soy protein antibodies
	Neogen	Veratox	Soy Flour Protein	2.5 (LOQ)	2.5	1 to 25 in PBS solution + additive	-
Walnut	Gen-Probe/Tepnel	Biokit	Walnut	0.25 (LOD)	2.4 - 120	1 to 10 in Tris (hydroxymethyl) aminomethane / sodium chloride / gelatin	Polyclonal antibody specifically detects walnut proteins

Table 2. Website Links for Detailed Information on the ELISA Commercial Kits Used for Determining Allergen Cross-Reactivity

Allergen	Company	Test kit	Detailed method (Website link)
Almonds	Gen-Probe/Tepnel	Biokit	http://www.tepnel.com/
	Neogen	Veratox	http://www.neogen.com/foodsafety/V_Product_List.asp?Test_Kit_Cat=203
	R-Biopharm	Ridascreen	http://www.r-biopharm.com/main.php?
Gluten	Gen-Probe/Tepnel	Biokit	http://www.tepnel.com/
	Neogen	Veratox	http://www.neogen.com/foodsafety/V_Product_List.asp?Test_Kit_Cat=203
	R-Biopharm	Ridascreen	http://www.r-biopharm.com/main.php?
Hazelnut	Gen-Probe/Tepnel	Biokit	http://www.tepnel.com/
	Neogen	Veratox	http://www.neogen.com/foodsafety/V_Product_List.asp?Test_Kit_Cat=203
	R-Biopharm	Ridascreen	http://www.r-biopharm.com/main.php?
Mustard	ELISA System	Mustard Seed P.	http://www.elisas.com.au/
	Neogen	Veratox	http://www.neogen.com/foodsafety/V_Product_List.asp?Test_Kit_Cat=203
	Sedum R&D	ELISA Kit	http://www.sedum-rd.cz/en/biotech/products/allergens-in-food/mustard-elisa-kit
Peanut	Gen-Probe/Tepnel	Biokit	http://www.tepnel.com/
	Neogen	Veratox	http://www.neogen.com/foodsafety/V_Product_List.asp?Test_Kit_Cat=203
	R-Biopharm	Ridascreen	http://www.r-biopharm.com/main.php?
Sesame	ELISA System	Sesame Residue	http://www.elisas.com.au/
	Gen-Probe/Tepnel	Biokit	http://www.tepnel.com/
	ELISA System	Soy Protein	http://www.elisas.com.au/
Soy	Neogen	Veratox	http://www.neogen.com/foodsafety/V_Product_List.asp?Test_Kit_Cat=203
Walnut	Gen-Probe/Tepnel	Biokit	http://www.tepnel.com/

**Table 3: Protein Hits Identified in the MASCOT Mass Spectrometric Database
with a MOWSE Score of 38 or Higher for the CDC Maria, C05041 and C05091**

Canary Seed Sample

Mass spectrometry protein hits (common name in bracket):

CDC Maria (*sample as is*)

AVE3_AVESA Avenin-3/Avena sativa (oats)
SSG1_AVESA 12S seed storage globulin 1/Avena sativa (oats)
AVEE_AVESA Avenin-E/Avena sativa (oats)
SSG2_AVESA 12S seed storage globulin 2/Avena sativa (oats)
GLUA2_ORYSJ Glutelin type-A 2/Oryza sativa subsp. japonica (rice)
G3PC_ANTMA Glyceraldehyde-3-phosphate dehydrogenase, cytosolic/Antirrhinum majus
(herbaceous perennial plant; common name – snapdragon)
ALF_CICAR Fructose-bisphosphate aldolase, cytoplasmic isozyme/Cicer arietinum (chickpea)
EF1A1_DAUCA Elongation factor 1-alpha/Daucus carota (wild carrot)
OLEO2_ORYSI Oleosin 18 kDa/Oryza sativa subsp. indica (rice)
PDI_WHEAT Protein disulfide-isomerase/Triticum aestivum (wheat)*
H4_ARATH Histone H4/Arabidopsis thaliana (annual flowering plant)
HSP70_MAIZE Heat shock 70 kDa protein/Zea mays (corn)

CDC Maria (*cleaned sample*)

AVE3_AVESA Avenin-3/Avena sativa (oats)
AVEE_AVESA Avenin-E/Avena sativa (oats)
SSG1_AVESA 12S seed storage globulin 1/Avena sativa (oats)
SSG2_AVESA 12S seed storage globulin 2/Avena sativa (oats)
GLUA2_ORYSJ Glutelin type-A 2/Oryza sativa subsp. japonica (rice)
H4_ARATH Histone H4/Arabidopsis thaliana (annual flowering plant)
H2B10_ORYSI Histone H2B.10/Oryza sativa subsp. indica (rice)
G3PE_MAIZE Glyceraldehyde-3-phosphate dehydrogenase, cytosolic 3/Zea mays (corn)
G3PX_HORVU Glyceraldehyde-3-phosphate dehydrogenase, cytosolic/Hordeum vulgare (barley)*
EM4_WHEAT Em protein H5/Triticum aestivum (wheat)*
SEEP_RAPSA Late seed maturation protein P8B6/Raphanus sativus (radish)

EMB1_DAUCA EMB-1 protein/Daucus carota (wild carrot)
 EMP1_ORYSJ Embryonic abundant protein 1/Oryza sativa subsp. japonica (rice)

C05041 (sample as is)

SSG1_AVESA 12S seed storage globulin 1/Avena sativa (oats)
 AVE3_AVESA Avenin-3/Avena sativa (oats)
 G3PE_MAIZE Glyceraldehyde-3-phosphate dehydrogenase, cytosolic 3/Zea mays (corn)
 GLUA2_ORYSJ Glutelin type-A 2/Oryza sativa subsp. japonica (rice)
 AVEE_AVESA Avenin-E/Avena sativa (oats)
 G3PX_HORVU Glyceraldehyde-3-phosphate dehydrogenase, cytosolic/Hordeum vulgare (barley)*
 HSP70_MAIZE Heat shock 70 kDa protein/Zea mays (corn)
 H2B10_ORYSI Histone H2B.10/Oryza sativa subsp. indica (rice)
 GLUB4_ORYSJ Glutelin type-B 4/Oryza sativa subsp. japonica (rice)
 GLUB1_ORYSJ Glutelin type-B 1/Oryza sativa subsp. japonica (rice)
 H4_ARATH Histone H4/Arabidopsis thaliana (annual flowering plant)
 ATPBM_HEVBR ATP synthase subunit beta, mitochondrial/Hevea brasiliensis (rubber tree)
 OLEO2_ORYSI Oleosin 18 kDa/Oryza sativa subsp. indica (rice)
 ATP9_MAIZE ATP synthase subunit 9, mitochondrial/Zea mays (corn)
 EF1A1_DAUCA Elongation factor 1-alpha/Daucus carota (wild carrot)
 ALF_CICAR Fructose-bisphosphate aldolase, cytoplasmic isozyme/Cicer arietinum (chickpea)

C05041 (cleaned sample)

AVE3_AVESA Avenin-3/Avena sativa (oats)
 AVEE_AVESA Avenin-E/Avena sativa (oats)
 SSG2_AVESA 12S seed storage globulin 2/Avena sativa (oats)
 H32_CICIN Histone H3.2/Cichorium intybus (chicory)
 H32_MEDSA Histone H3.2/Medicago sativa (alfafa)
 GLUA2_ORYSJ Glutelin type-A 2/Oryza sativa subsp. japonica (rice)
 G3PE_MAIZE Glyceraldehyde-3-phosphate dehydrogenase, cytosolic 3/Zea mays (corn)
 UP01_VITRO Unknown protein 1 (Fragment)/Vitis rotundifolia (grapevine species)
 G3PX_HORVU Glyceraldehyde-3-phosphate dehydrogenase, cytosolic/Hordeum vulgare (barley)*
 EF2_BETVU Elongation factor 2/Beta vulgaris (beet)
 H2B10_ORYSI Histone H2B.10/Oryza sativa subsp. indica (rice)
 SEEP_RAPSA Late seed maturation protein P8B6/Raphanus sativus (radish)
 EMB1_DAUCA EMB-1 protein/Daucus carota (wild carrot)
 H4_ARATH Histone H4/Arabidopsis thaliana (annual flowering plant)
 SUSY_ALNGL Sucrose synthase/Alnus glutinosa (alder)
 SUS1_ORYSJ Sucrose synthase 1/Oryza sativa subsp. japonica (rice)

GRDH_DAUCA Glucose and ribitol dehydrogenase/Daucus carota (wild carrot)
 GRDH_ORYSJ Glucose and ribitol dehydrogenase homolog/Oryza sativa subsp. japonica (rice)
 EMP1_ORYSJ Embryonic abundant protein 1/Oryza sativa subsp. japonica (rice)
 PROF2_SOLLC Profilin-2/Solanum lycopersicum (tomato)
 REHYA_ORYSI 1-Cys peroxiredoxin A/Oryza sativa subsp. indica (rice)
 EF1A1_DAUCA Elongation factor 1-alpha/Daucus carota (wild carrot)

C05091 (sample as is)

SSG1_AVESA 12S seed storage globulin 1/Avena sativa (oats)
 AVE3_AVESA Avenin-3/Avena sativa (oats)
 SSG2_AVESA 12S seed storage globulin 2/Avena sativa (oats)
 AVEE_AVESA Avenin-E/Avena sativa (oats)
 GLUA2_ORYSJ Glutelin type-A 2/Oryza sativa subsp. japonica (rice)
 G3PC_ANTMA Glyceraldehyde-3-phosphate dehydrogenase, cytosolic/Antirrhinum majus
 (herbaceous perennial plant; common name – snapdragon)
 EF1A1_DAUCA Elongation factor 1-alpha/Daucus carota (wild carrot)

C05091 (cleaned sample)

AVE3_AVESA Avenin-3/Avena sativa (oats)
 SSG1_AVESA 12S seed storage globulin 1/Avena sativa (oats)
 SSG2_AVESA 12S seed storage globulin 2/Avena sativa (oats)
 AVEE_AVESA Avenin-E/Avena sativa (oats)
 GLUA2_ORYSJ Glutelin type-A 2/Oryza sativa subsp. japonica (rice)
 H4_ARATH Histone H4/Arabidopsis thaliana (annual flowering plant)
 EF1A1_DAUCA Elongation factor 1-alpha/Daucus carota (wild carrot)
 OLEO2_ORYSI Oleosin 18 kDa/Oryza sativa subsp. indica (rice)
 HSP70_MAIZE Heat shock 70 kDa protein/Zea mays (corn)
 SUS1_ORYSJ Sucrose synthase 1/Oryza sativa subsp. japonica (rice)
 SEEP_RAPSA Late seed maturation protein P8B6/Raphanus sativus (radish)
 EMB1_DAUCA EMB-1 protein/Daucus carota (wild carrot)

(* wheat and/or barley)

Table 4. List of Wheat, Non- Gluten Cereals and Canary Seed Proteins Identified by LC/MS Analysis of Excised Bands from 1-DE, Selected by IgE Abs in Immunoblotting

Band No. ^a	Accession numbers	Protein identified	Sequence coverage ^b (%)	Mascot Score ^c	Theoretical mass (kDa)
W1	B7U6L4	Wheat globulin 3	24	50	66
	A9YSK3	HMW glutenin subunit	14	52	70
W2	B2Y2R5	LMW glutenin subunit	12	52	45
W3	P17314	Wheat alpha amylase/ trypsin inhibitor	20	52	18
W4	P01085 (Q19A44)	Wheat alpha amylase/ trypsin inhibitor	95	52	13
W5	Q5UHH8	Wheat alpha amylase/ trypsin inhibitor	90	52	13
O1	Q38780	AVESA 11S globulin (Avena Sativa=Oat)	38	53	59
O2	P14812	AVESA 12S globulin 2	33	53	59

O3	O49258	AVESA 12S globulin	25	53	58
O4	O49257	AVESA 12S globulin	14	51	53
O5	Q38780	AVESA 11S globulin	28	53	59
Q1-Q5	Q06AW2	CHEQI 11S globulin A/B	13	52	54
	(Q06AW1)				
M1-M2	DX3508	Granule bound starch	72	52	58
	DX3516	synthase I. (GBSSI)			
T1-T2		Not identified by Mascot			
S1	A5Y380	SORBI-GBSSI (from Sorghum bicolor)	55	54	66
B1	O23878	FAGES 13S globulin	13	52	65
		(Fagopyrum esculentum)			
B2	O23878	FAGES 13S globulin	1 7.5	52	65
		(Fagopyrum esculentum)			
B3	A5H1X6	FAGES Vicilin-like protein	24	51	16
	Q6QJL1	FAGES Vicilin-like protein	17	52	16
C1	ND	<u>Significant homology with:</u>			
		Wheat (GBSS1=granule-bound starch synthase 1)	13	63.7	68

		Maize	10	105.0	73
		Rice	7.7	80.5	61
		Oat	6.9	45.5	62
		Sorghum	5.3	62.2	55
		Barley	2.5	53.3	72
C2	ND	<u>Significant homology with:</u>			
		<u>Rice</u>	20	49.9	66
		<u>Oat</u>	13	45.5	53
		<u>Barley</u>	7.6	53.3	61
		<u>Maize</u>	4.2	53.0	68
		<u>Wheat (GBSS1)</u>	2.8	54.1	64
C3	ND	<u>Significant homology with:</u>			
		<u>Barley</u>	29	73.3	59
		<u>Maize</u>	14	70.3	52
		<u>Oat</u>	11	52.1	59
		<u>Wheat (GBSS1)</u>	7.6	52.6	57
		<u>Rice</u>	6.5	84.2	57
C4	ND	<u>Significant homology with:</u>			
		<u>Rice</u>	16	57.0	45

		<u>Maize</u>	11	70.0	43
		<u>Oat</u>	10	75.0	58
		<u>Wheat (GBSS1)</u>	6.8	78.6	57
		<u>Sorghum</u>	5.7	58.2	46
C5	ND	<u>Significant homology with:</u>			
		<u>Oat</u>	12	46	39
		<u>Barley</u>	12	58.4	33
C6	ND	<u>Significant homology with:</u>			
		<u>Barley</u>	17	59.6	27
		<u>Oat</u>	16	39.7	25
		<u>Rice</u>	16	57.0	17
C7	ND	<u>Significant homology with:</u>			
		<u>Rye</u>	29	82.1	10
		<u>Maize</u>	17	74	17
		<u>Oat</u>	14	51.5	15

^aBand numbers correspond to the gel in Figure 5; ^bSequence coverage b = Total percentage of proteins amino acid sequence covered by peptides in MS/MS analysis. ^cMASCOT score c = >40 indicate identification or extensive homology; HMW - high molecular weight; LMW - low molecular weight. ND: not determined (there are no canary seed protein sequences deposited in the NCBI dbEST).

**Table 5. ELISA Results of Glabrous Canary Seed Analysed as is (Uncleaned)
and Using Quantity Indicated in Kit Instructions**

Allergen	Company	Test kit	ELISA results		
			CDC Maria	C05041	C05091
			ppm	ppm	ppm
Almond	Gen-Probe/Tepnel	Biokit	< LOD	< LOD	< LOD
	Neogen	Veratox	< LOQ	< LOQ	< LOQ
	R-Biopharm	Ridascreen	< LOD	< LOD	< LOD
Gluten	Gen-Probe/Tepnel	Biokit	< LOD	< LOD	< LOD
	Neogen	Veratox	5.4 ($\cong$ LOQ)	< LOQ	< LOQ
	R-Biopharm	Ridascreen	< LOD	< LOD	< LOD
Hazelnut	Gen-Probe/Tepnel	Biokit	< LOD	< LOD	< LOD
	Neogen	Veratox	< LOQ	< LOQ	< LOQ
	R-Biopharm	Ridascreen	< LOD	< LOD	< LOD
Mustard	ELISA System	ELISA	13.0 a, A	0.7 b	6.9 c
	Sedium R&D	ELISA	< LOD	< LOD	< LOD
	Neogen	Veratox	36.5 B	< LOQ	< LOQ
Peanut	Gen-Probe/Tepnel	Biokit	< LOD	< LOD	< LOD
	Neogen	Veratox	< LOQ	< LOQ	< LOQ
	R-Biopharm	Ridascreen	< LOD	< LOD	< LOD
Sesame	ELISA System	ELISA	< LOQ	< LOQ	< LOQ
	Gen-Probe/Tepnel	Biokit	< LOD	< LOD	< LOD
Soy	ELISA System	ELISA	< LOQ	< LOQ	< LOQ
	Neogen	Veratox	< LOQ	< LOQ	< LOQ
Walnut	Gen-Probe/Tepnel	Biokit	< LOD	< LOD	< LOD

LOD: Limit of detection; LOQ: Limit of quantification.

Small letters (a, b, c) within a row indicate significant ($p < 0.05$) differences within the canary seed composites.

Capital letters (A and B) within a column indicate significant ($p < 0.05$) differences between the ELISA kits for a particular variety.

**Table 6. ELISA Results of Glabrous Canary seed Analysed as is (Uncleaned)
and Using Three Times the Quantity Indicated in Kit Instructions**

Allergen	Company	Test kit	ELISA results		
			CDC Maria	C05041	C05091
Almond	Gen-Probe/Tepnel	Biokit	< LOD	< LOD	< LOD
	Neogen	Veratox	< LOQ	< LOQ	< LOQ
	R-Biopharm	Ridascreen	< LOD	< LOD	< LOD
Gluten	Gen-Probe/Tepnel	Biokit	< LOD	< LOD	< LOD
	Neogen	Veratox	< LOQ	< LOQ	< LOQ
	R-Biopharm	Ridascreen	< LOD	< LOD	< LOD
Hazelnut	Gen-Probe/Tepnel	Biokit	< LOD	< LOD	< LOD
	Neogen	Veratox	< LOQ	< LOQ	< LOQ
	R-Biopharm	Ridascreen	< LOD	< LOD	< LOD
Mustard	ELISA System	ELISA	< LOQ	< LOQ	< LOQ
	Sedium R&D	ELISA	< LOD	< LOD	< LOD
	Neogen	Veratox	< LOQ	< LOQ	< LOQ
Peanut	Gen-Probe/Tepnel	Biokit	< LOD	< LOD	< LOD
	Neogen	Veratox	< LOQ	< LOQ	< LOQ
	R-Biopharm	Ridascreen	< LOD	< LOD	< LOD
Sesame	ELISA System	ELISA	< LOQ	< LOQ	< LOQ
	Gen-Probe/Tepnel	Biokit	< LOD	< LOD	< LOD
Soy	ELISA System	ELISA	< LOQ	< LOQ	< LOQ
	Neogen	Veratox	< LOQ	< LOQ	< LOQ
Walnut	Gen-Probe/Tepnel	Biokit	< LOD	< LOD	< LOD

LOD: Limit of detection; LOQ: Limit of quantification.

Table 7. ELISA Results for Glabrous Canary seed Analysed after Cleaning and Using One and Three Times the Quantity Indicated in Kit Instructions

Allergen	Company	Test kit	ELISA results		
			CDC Maria	C05041	C05091
A. Canary seed analysed after cleaning and using quantity indicated in kit instructions					
Gluten	Neogen	Veratox	< LOQ	< LOQ	< LOQ
		Veratox R5	< LOQ	< LOQ	< LOQ
Mustard	ELISA System	ELISA	< LOQ	< LOQ	< LOQ
	Neogen	Veratox	< LOQ	< LOQ	< LOQ
B. Canary seed analysed after cleaning and using three times the quantity indicated in kit instructions					
Gluten	Neogen	Veratox	< LOQ	< LOQ	< LOQ
Mustard	ELISA System	ELISA	< LOQ	< LOQ	< LOQ
	Neogen	Veratox	< LOQ	< LOQ	< LOQ

LOQ: Limit of quantification

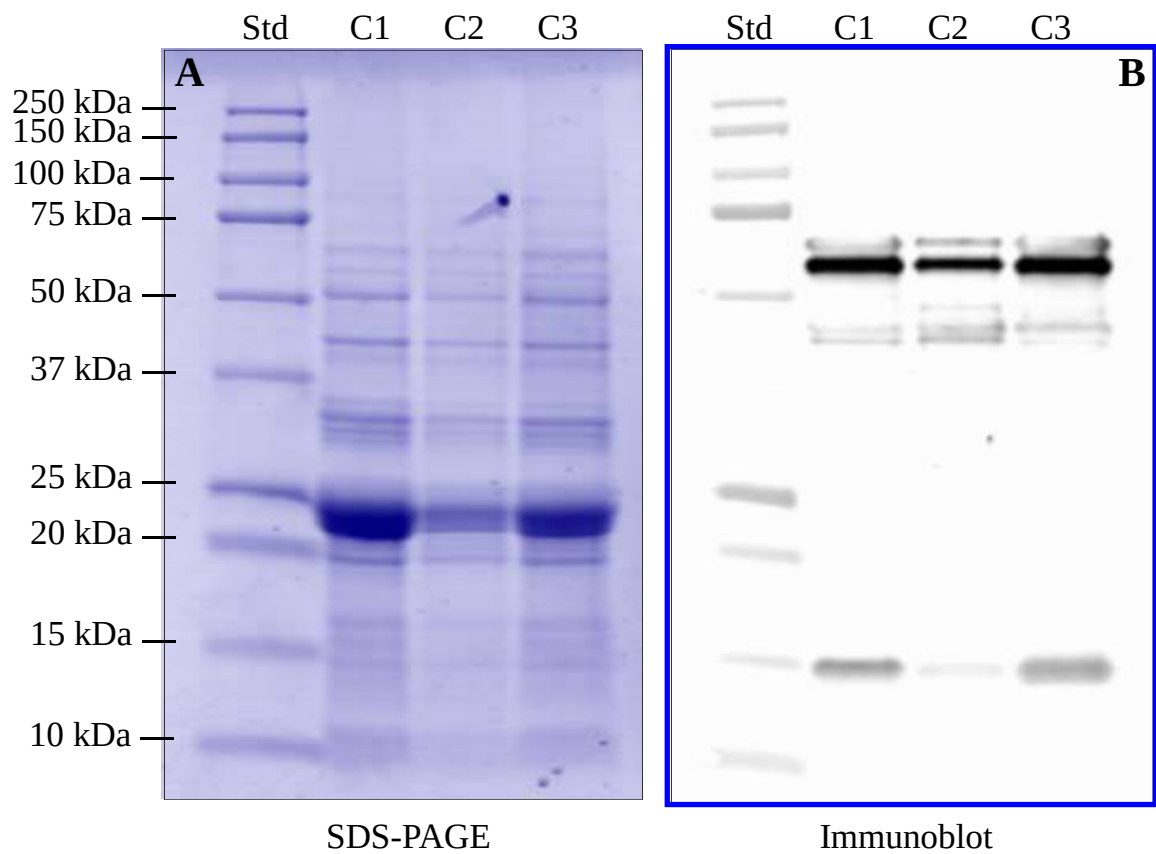
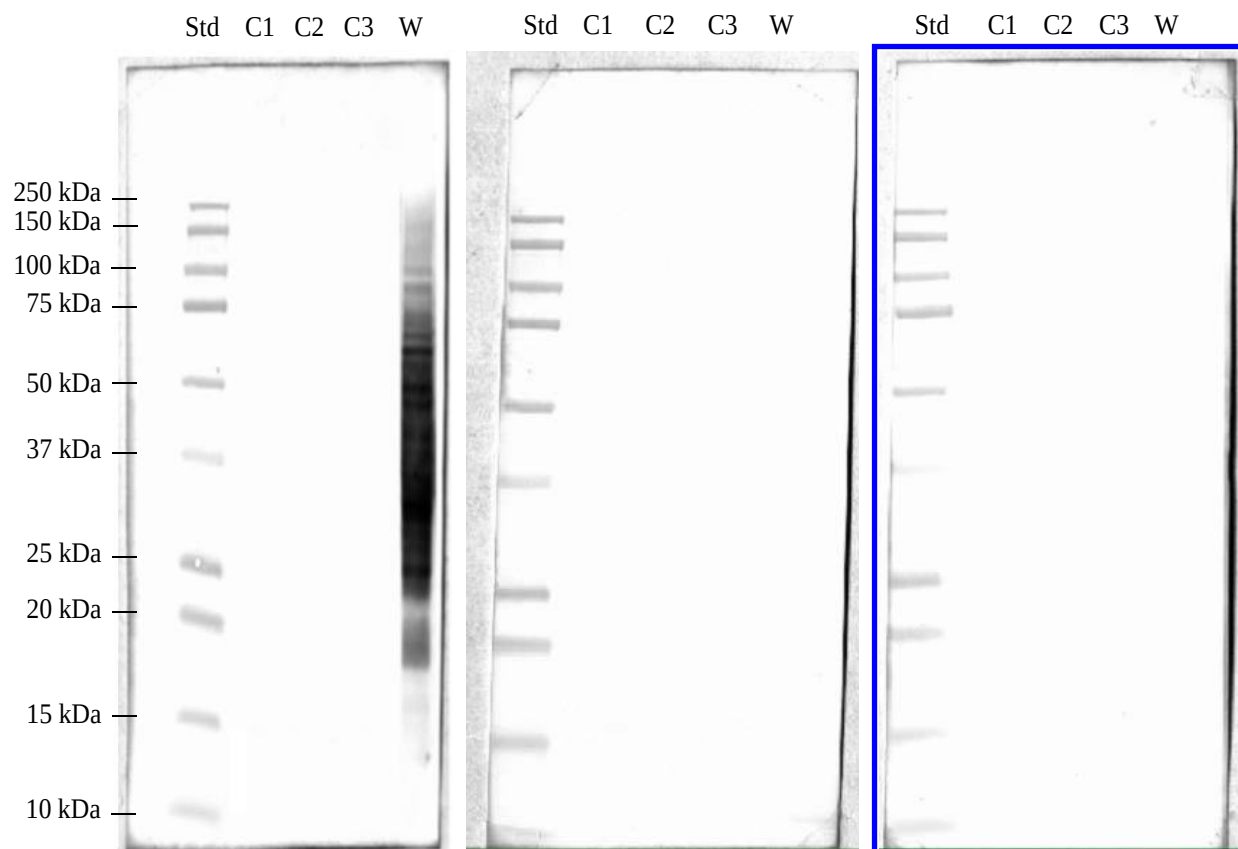


Figure 1

**Figure 2**


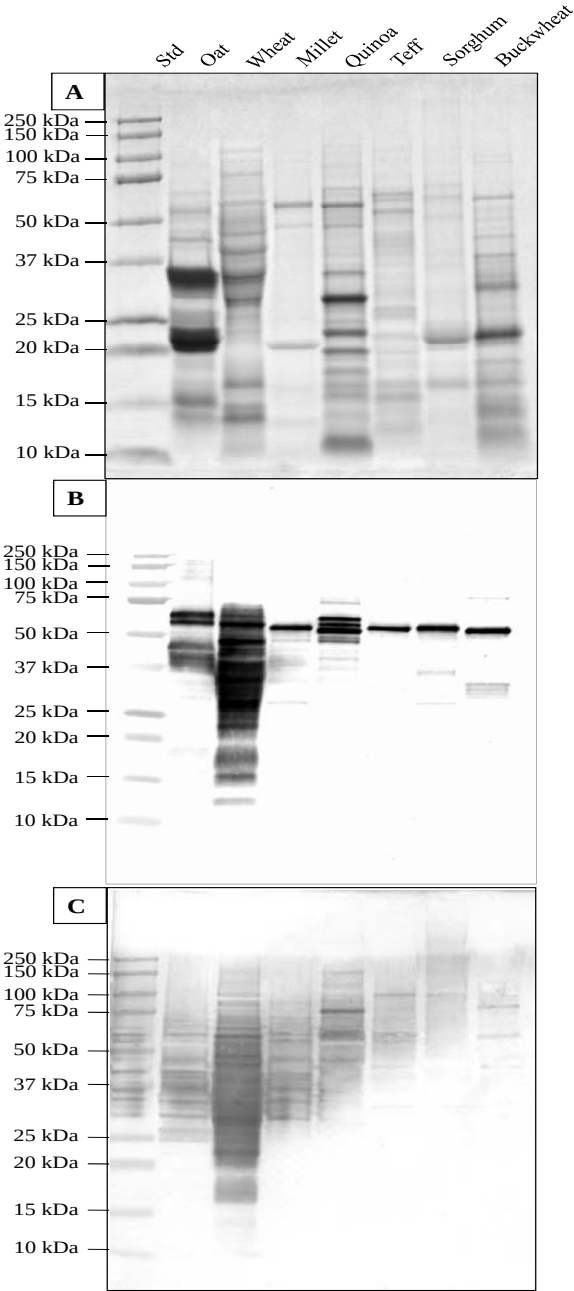


Figure 3

**Figure 4**

TOC GRAPHIC

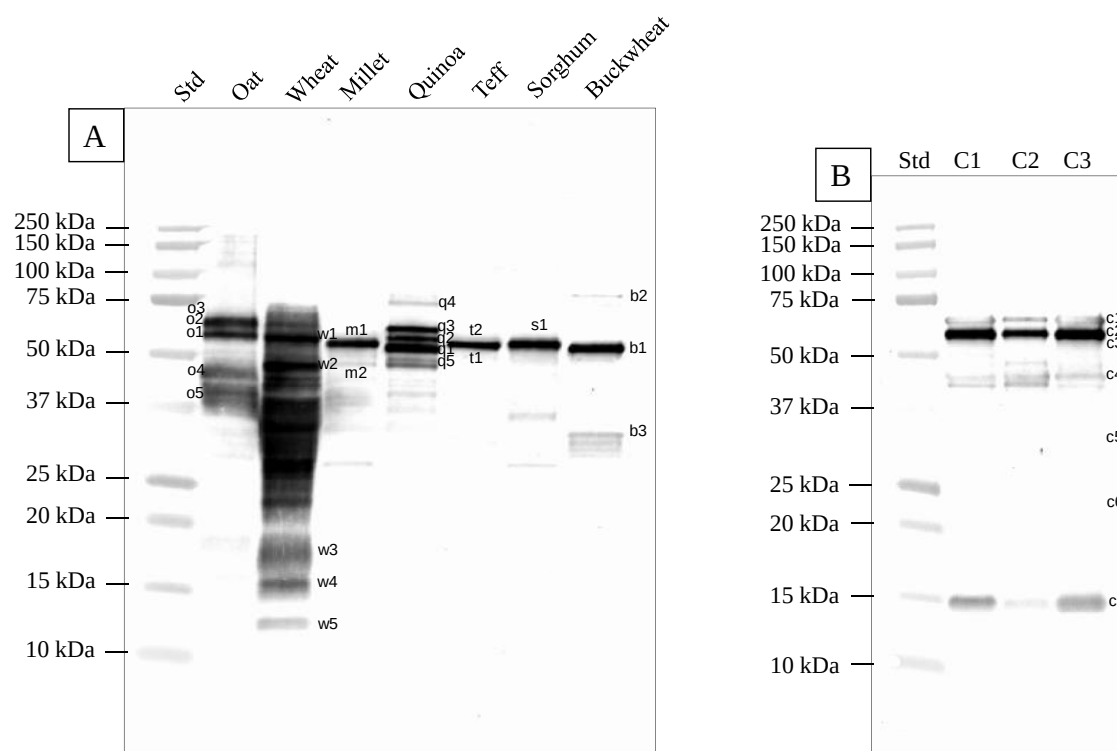